

Product: Rabbit Anti-Histone H3, affinity purified, react with histone H3 from various species including human, mouse, rat, bovine etc.

Catalog Number: AC0003

Lot Number: 10923

Size: 0.1ml affinity purified Abs, 1mg/ml in PBS with 0.02% Na Azide


Immunogen: Histone H3 peptide C-terminus (MPKDIQLARRIGERA).

Specification: AC0003 binds to various histone H3 proteins of many different species including human, mouse, rat, bovine, rabbit, xenopus etc in Western Blot analysis.

Suggested Usage: The product can be used at 1:1,000 to 1:5,000 for immunoblot testing. The product can be used in ELISA, Western Blot and other immuno-assays. User needs to optimize the concentrations in assay. The product can be used in combination with GeneTel's anti-IgY conjugates such as RCYHRP (anti-IgY HRP) and RCYAP (anti-IgY AP) Western Blot analysis.

Short Term Storage: 2-8°C for one week.

Long Term Storage: -20°C as an undiluted liquid for up to two year.
Please feel free to call us if you have any question regarding to this product or any other GeneTel's product.

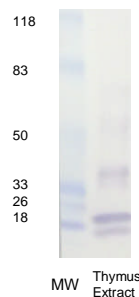

Western Blot analysis of rabbit anti-histone (AC0003) shows that the Abs reacts with various H3 proteins on calf thymus protein extract.