

**GeneTel Laboratories, LLC** *Leader in Proteomics Abs and Abs Arrays*  
1202 Ann Street, Madison, WI 53713 Phone 877-248-4316 Fax 608-257-1385 Email info@genetel-lab.com

**Product:** Chicken Anti-human BAG family molecular chaperone regulator 1 (BAG1), total IgY.

**Catalog Number:** AP0007

**Lot Number:** 70502

**Size:** 0.5ml Abs, 1mg/ml in PBS with 0.02% Na Azide

**Immunogen:** synthetic peptide


**Specification:** AP0007 reveals a band at the expected position (31KD) in a Western Blot analysis.

**Suggested Usage:** The product can be used at 1:500 to 1:5000 for immunoblot testing. User needs to optimize the concentrations in assay. The product can be used in combination with GeneTel's anti-chicken IgY conjugates in Western Blot analysis. This product may react with homologues from other mammalian species.

**Short Term Storage:** 2-8°C for one week.

**Long Term Storage:** -20°C as an undiluted liquid for up to two years.

Please feel free to call us if you have any question regarding to this product or any other GeneTel's products.

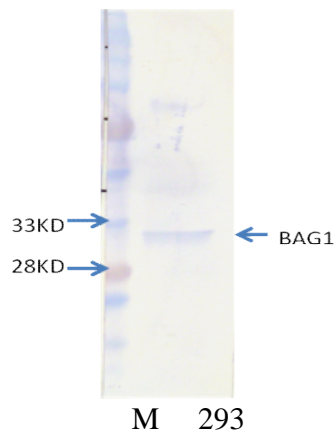

**Western Blot analysis of AP0007**

M: protein markers

Lanes 293: 293 lysate

1<sup>st</sup> antibody AP0007 (1:2000)

2<sup>nd</sup> antibody: Rabbit anti-chicken IgY-AP

Developed with chromogenic substrate.