

GeneTel Laboratories, LLC *Leader in Proteomics Abs and Abs Arrays*
1202 Ann Street, Madison, WI 53713 Phone 877-248-4316 Fax 608-257-1385 Email info@genetel-lab.com

Product: Chicken Anti-human caspase 3, total IgY
Catalog Number: AP0036

Lot Number: 71802

Size: 0.5ml Abs, 1mg/ml in PBS with 0.02% Na Azide

Immunogen: synthetic peptide


Specification: AP0036 reveals a band at the expected position in a Western Blot analysis.

Suggested Usage: The product can be used at 1:1000 to 1:5000 for immunoblot testing. User needs to optimize the concentrations in assay. The product can be used in combination with GeneTel's anti-chicken IgY conjugates in Western Blot analysis. This product may react with homologues from other mammalian species.

Short Term Storage: 2-8°C for one week.

Long Term Storage: -20°C as an undiluted liquid for up to two years.

Please feel free to call us if you have any question regarding to this product or any other GeneTel's products.

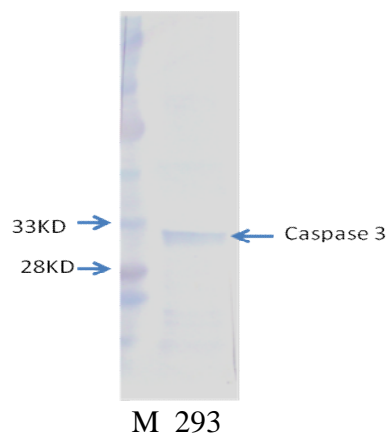

Western Blot analysis of AP0036

M: protein markers

Lanes 293: 293cell lysate

1st antibody AP0036 (1:5000)

2nd antibody: Rabbit anti-chicken IgY-AP

Developed with chromogenic substrate.