

Genetel Laboratories, LLC

Leader in Proteomics Abs and Abs Arrays

1202 Ann Street, Madison, WI 53713 Phone 877-248-4316 Fax 608-257-1385 Email info@genetel-lab.com

Product: Rabbit anti-GFAP, affinity purified pAbs against glial fibrillary acid protein (GFAP), an astrocyte marker of central nervous system. Human, mouse and rat specific.

Catalog Number: CN0021

Lot Number: 120530

Size: 50ug in 50ul PBS

Immunogen: GFAP fusion protein


Suggested Usage: The product can be used at 1:1,000 to 1:2,000 for immunoblot testing or 1:500-1:1000 for IHC or ICC. The product can be used in combination with GeneTel's conjugates such as DRAP, DRCY3, DRFITC or GRHRP.

Short Term Storage: 2-8°C for one week.

Long Term Storage: -20°C as an undiluted liquid for up to two year.

Note: For in vitro research use only, not for diagnostic or therapeutic use.

Please feel free to call us if you have any question regarding to this product or any other GeneTel's product.

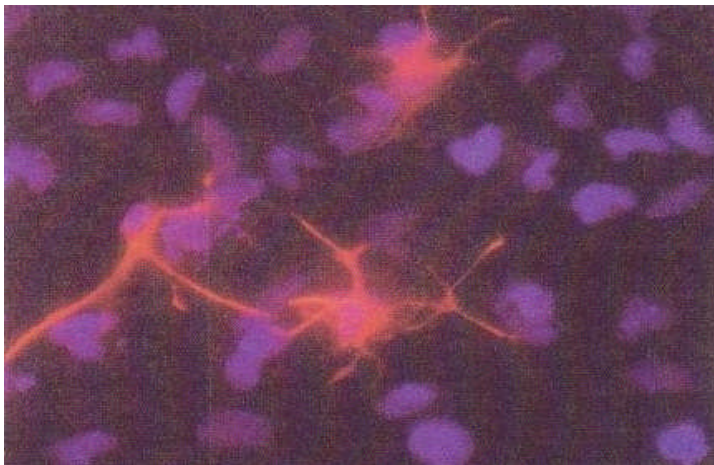

Anti-GFAP Immunostaining of astrocytes in rat neural progenitor cell culture. CN0021 with anti-rabbit IgG Cy3 (DRCY3, red) and DAPI (blue) were used.