

Product: Chicken Anti-GST protein, affinity purified.

Catalog Number: GST

Lot Number: 71112

Size: 0.1ml, 0.5mg/ml in PBS with 0.02% Na Azide

Immunogen: Purified GST protein.


Specification: The antibodies bind specifically to the GST protein in cell lysate containing GST or GST fusion protein in a western blot analysis and other immunoassays.

Suggested Usage: The product can be used at 1:2,000 to 1:5,000 for both immunoblot and ELISA testing. User needs to try out the titrations for optimization in immunocytochemistry (ICC) and immuno-histochemistry (IHC) assay. The product can be used in combination with GeneTel's anti-IgY conjugates such as RCYHRP (anti-IgY HRP) and RCYAP (anti-IgY AP) for ELISA and Western Blot analysis, and RCYRP or RCYFITC for ICC/IHC assay.

Short Term Storage: 2-8°C for one week.

Long Term Storage: -20°C as an undiluted liquid for up to two year.

Please feel free to call us if you have any question regarding to this product or any other GeneTel's product.

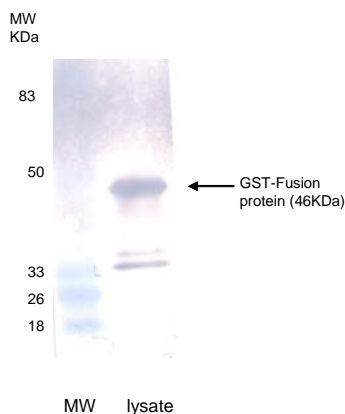

Western Blot analysis of cell lysate containing a GST fusion protein with the chicken anti-GST at a 1:4,000 dilution.