

Product: Chicken Anti-DHHC5, zinc finger DHHC domain containing 5, affinity purified. Human, mouse and rat specific.

Catalog Number: OA0019

Lot Number: 80229

Size: 0.1mg in 0.1ml PBS with 0.02% Na Azide


Immunogen: Human and mouse DHHC5 peptide aa245-261.

Specification: OA0019 binds to a strong band at about 80KDa on HEK293 cell line, the apparent MW may be bigger (about 110KDa) in liver protein lysate. The product may cross react with DHHC5 from other mammalian species.

Suggested Usage: The product can be used at 1:2,000 to 1:5,000 for immunoblot testing and 1:50-1:200 for immunoprecipitation.. User needs to optimize the concentrations in assay. The product can be used in combination with GeneTel's anti-IgY conjugates such as RCYHRP (anti-IgY HRP) and RCYAP (anti-IgY AP) Western Blot analysis or BY0001 (bovine anti-IgY resin) for immunoprecipitation.

Short Term Storage: 2-8°C for one week.

Long Term Storage: -20°C as an undiluted liquid for up to two year.
Please feel free to call us if you have any question regarding to this product or any other GeneTel's product.

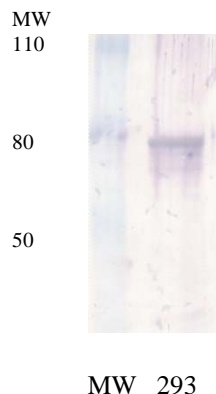

Western Blot analysis of OA0019 chicken anti-DHHC5 (1:5000) on HEK293 cell protein lysate.