

Product: Rabbit Anti-UBC12, Ubiquitin Conjugating Enzyme E2M , affinity purified. Human, mouse and rat specific.

Catalog Number: OA0025

Lot Number: 80821

Size: 0.1ml affinity purified Abs, 0.45mg/ml in PBS with 0.02% Na Azide


Immunogen: Human UBC12 C-terminus peptide.

Specification: OA0025 binds to a protein band at a MW of about 22KDa (UBC12) in human or mouse liver and spleen protein lysate in a Western Blot analysis. The product may cross react with UBC12 from other mammalian species.

Suggested Usage: The product can be used at 1:1,000 to 1:5,000 for immunoblot testing. User needs to optimize the concentrations in assay. The product can be used in combination with GeneTel's anti-rabbit IgG conjugates such as GRHRP in Western Blot analysis.

Short Term Storage: 2-8°C for one week.

Long Term Storage: -20°C as an undiluted liquid for up to two year. Please feel free to call us if you have any question regarding to this product or any other GeneTel's product.

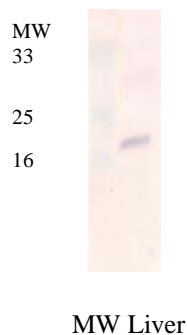

Western Blot analysis of OA0025 rabbit anti-UBC12 (1:1000) on mouse liver protein lysate (40ug). The Abs binds to UBC12, a protein at about 22 KDa.