

Product: Rabbit Anti-Sirt3 (hSIRT3), Silent Information Regulator 3 or NAD-dependent deacetylase sirtuin 3, affinity purified, human specific

Catalog Number: OA0076

Lot Number: 10802

Size: 0.1mg affinity purified Abs at 1mg/ml in PBS with 0.02% Na Azide

Immunogen: Human hSIRT3 C-terminus peptide (EEMRDLVQRETGKLDGPKD).


Specification: OA0076 binds to a protein band at a MW of about 28KDa of human cell mitochondrial protein preparation in a Western Blot analysis.

Suggested Usage: The product can be used at 1:2,000 to 1:5,000 for immunoblot testing. User needs to optimize the concentrations in assay. The product can be used in combination with GeneTel's anti-IgY conjugates such as RCYHRP (anti-IgY HRP) and RCYAP (anti-IgY AP) Western Blot analysis.

Short Term Storage: 2-8°C for one week.

Long Term Storage: -20°C as an undiluted liquid for up to two year.

Please feel free to call us if you have any question regarding to this product or any other GeneTel's product.

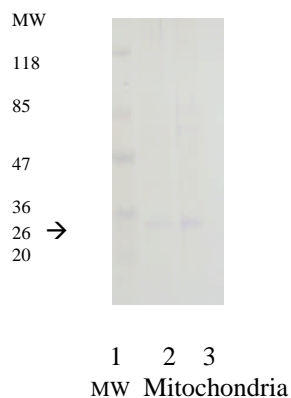

Western Blot analysis of OA0076 rabbit anti-hSIRT3 (1:5,000) on two human cell mitochondrial protein preparations (lane 2&3, 30ug protein). The Abs binds to the human Sirt3, a protein at about 28 KDa molecular weight.