

**GeneTel Laboratories, LLC** *Leader in Proteomics Abs and Abs Arrays*  
1202 Ann Street, Madison, WI 53713 Phone 877-248-4316 Fax 608-257-1385 Email info@genetel-lab.com

**Product:** Chicken Anti-human PDE9A, cGMP-specific 3',5'-cyclic phosphodiesterase type 9, total IgY

**Catalog Number:** OA0085

**Lot Number:** 120910

**Size:** 0.2ml Abs, 1mg/ml in PBS with 0.02% Na Azide

**Immunogen:** recombinant protein


**Specification:** OA0085 reveals a specific band with a MW of ~60KD in a Western Blot analysis, as expected.

**Suggested Usage:** The product can be used at 1:200 to 1:1000 for immunoblot testing. User needs to optimize the concentrations in assay. The product can be used in combination with GeneTel's anti-chicken IgY conjugates in Western Blot analysis. This product may react with homologues from other mammalian species.

**Short Term Storage:** 2-8°C for one week.

**Long Term Storage:** -20°C as an undiluted liquid for up to two years.

Please feel free to call us if you have any question regarding to this product or any other GeneTel's products.

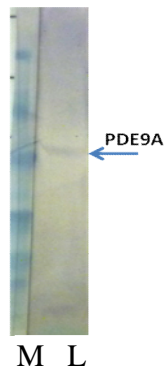

### Western Blot analysis of OA0085

M: protein markers

L: human lung tissue cell free extract

1<sup>st</sup> antibody: OA0085 (5ug/ml)

2<sup>nd</sup> antibody: Rabbit anti-chicken IgY-AP

Developed with chromogenic substrate.