

GeneTel Laboratories, LLC *Leader in Proteomics Abs and Abs Arrays*
1202 Ann Street, Madison, WI 53713 Phone 877-248-4316 Fax 608-257-1385 Email info@genetel-lab.com

Product: Chicken Anti-human chromosome 21 orf97.1, total IgY
Catalog Number: OA0088

Lot Number: 120910

Size: 0.2ml Abs, 1mg/ml in PBS with 0.02% Na Azide

Immunogen: recombinant protein


Specification: OA0088 reveals a specific band with a MW of about 30 KD in a Western Blot analysis, as predicted.

Suggested Usage: The product can be used at 1:200 to 1:1000 for immunoblot testing. User needs to optimize the concentrations in assay. The product can be used in combination with GeneTel's anti-chicken IgY conjugates in Western Blot analysis. This product may react with homologues from other mammalian species.

Short Term Storage: 2-8°C for one week.

Long Term Storage: -20°C as an undiluted liquid for up to two years.

Please feel free to call us if you have any question regarding to this product or any other GeneTel's products.

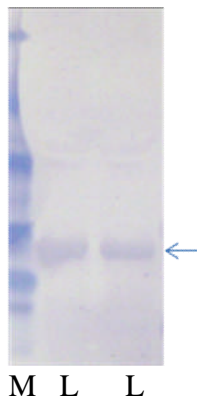

Western Blot analysis of OA0088

M: protein markers

Lanes L: human lung tissue homogenate (0.2mg total protein loading)

1st antibody: OA0088 (5ug/ml)

2nd antibody: Rabbit anti-chicken IgY-AP

Developed with chromogenic substrate.