

GeneTel Laboratories, LLC *Leader in Proteomics Abs and Abs Arrays*
1202 Ann Street, Madison, WI 53713 Phone 877-248-4316 Fax 608-257-1385 Email info@genetel-lab.com

Product: Rabbit polyclonal antibody against human ICAM-1, antigen affinity purified

Cat. No. OA0105

Lot Number: 2

Concentration: 0.1mg/0.1ml in PBS/Trehalose


Immunogen: recombinant human ICAM-1 protein

Specification: In western blot, this antibody specifically recognizes recombinant human ICAM-1 protein

Suggested Usage: The product can be used for ELISA and immuno blotting. Dilution of 1:1000-1:5000 is suggested for Western blot, the optimal concentration, however, needs to be experimentally determined by users.

Short Term Storage: 2-8°C for one week.

Long Term Storage: -20°C as an undiluted liquid for up to two years.

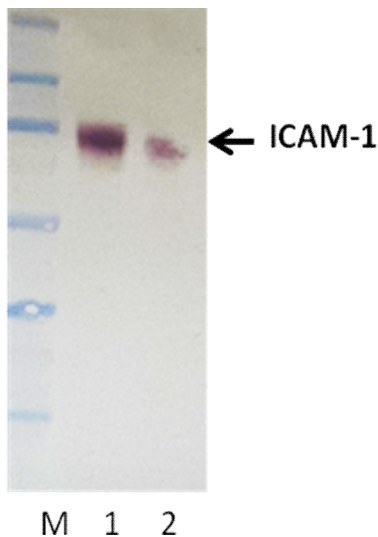

Western Blot analysis on OA0105

M: Molecular weight marker

1: 50ng rhIP10 protein

2: 50ng rhIP10 protein spiked in 293 cell lysate (25ug protein)

1st antibody used at 0.2ug/ml

2nd antibody used at 1:5000 (AP conjugated donkey anti-rabbit IgG)

Color was developed with chromogenic substrate.

* Two variants of recombinant IP10 expressed and purified from E. coli.