



GeneTel Laboratories, LLC

Leader in Proteomics Abs and Abs Arrays

1202 Ann Street, Madison, WI 53713 Phone 877-248-4316 Fax 608-257-1385 Email info@genetel-lab.com

Product: Chicken Anti-PPAR γ (Peroxisome Proliferative Activated Receptor gamma), affinity purified.

Catalog Number: OB0010

Lot Number: 61024

Size: 0.1ml PBS with 0.02% Na Azide, 0.5mg/ml IgY Antibodies.

Immunogen: Human PPAR γ peptide (YYSEKTQLYNKPHEE).


Specification: The antibody reacts specifically with PPAR γ protein of human and may cross-react with PPAR γ protein of other mammalian species such as swine, mouse and rat.

Suggested Usage: The product can be used at 1:2000-1:5000 dilutions for Western Blot analysis and for other immuno-assays. Customers need to optimize the dilution when using the Abs in IHC or ICC.

Short Term Storage: 2-8°C for one week.

Long Term Storage: -20°C as an undiluted liquid for up to two year.

Note: For in vitro research use only, not for diagnostic or therapeutic use.

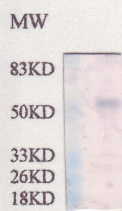

Western Blot analysis of OB0010 Chicken Anti-PPAR γ Abs (1:2000) on human adipose protein lysate (50ug).