



GeneTel Laboratories, LLC

Leader in Proteomics Abs and Abs Arrays

1202 Ann Street, Madison, WI 53713 Phone 877-248-4316 Fax 608-257-1385 Email info@genetel-lab.com

Product: Chicken Anti-Human ALBP (Adipocyte Lipid or Fatty Acid Binding Protein, or aP2), affinity purified.

Catalog Number: OB0011

Lot Number: 61102

Size: 0.1ml PBS with 0.02% Na Azide, 0.5mg/ml IgY Antibodies.

Immunogen: Human ALBP peptide (KWDGKSTTIKRKRED).


Specification: The antibody reacts specifically with human ALBP (about 16KD) and may cross react with ALBP of other mammalian species.

Suggested Usage: The product can be used at 1:2000-1:5000 dilutions for Western Blot analysis and for other immuno-assays. Customers need to optimize the dilution when using the Abs in IHC or ICC.

Short Term Storage: 2-8°C for one week.

Long Term Storage: -20°C as an undiluted liquid for up to two year.

Note: For in vitro research use only, not for diagnostic or therapeutic use.

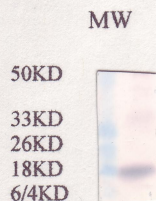

Western Blot analysis of OB0011 Chicken Anti- ALBP Abs (1:2000) on human adipose tissue protein lysate (50ug).