



**GeneTel Laboratories, LLC**

*Leader in Proteomics Abs and Abs Arrays*

1202 Ann Street, Madison, WI 53713 Phone 877-248-4316 Fax 608-257-1385 Email [info@genetel-lab.com](mailto:info@genetel-lab.com)

**Product:** Chicken Anti-Mouse HSL (Hormone Sensitive Lipase, or LIPE), affinity purified.

**Catalog Number:** OB0016

**Lot Number:** 70116

**Size:** 0.1ml PBS with 0.02% Na Azide, 0.5mg/ml IgY Antibodies.

**Immunogen:** Mouse HSL peptide (LSGRKPQKTTSPТАES).


**Specification:** The antibody reacts specifically with mouse HSL (about 83KD) of mouse adipose protein lysate in western blot analysis, and may cross-react with an additional weak band at about 40KD.

**Suggested Usage:** The product can be used at 1:2000-1:5000 dilutions for Western Blot analysis and for other immuno-assays. Customers need to optimize the dilution when using the Abs in IHC or ICC.

**Short Term Storage:** 2-8°C for one week.

**Long Term Storage:** -20°C as an undiluted liquid for up to two year.

**Note:** For in vitro research use only, not for diagnostic or therapeutic use.

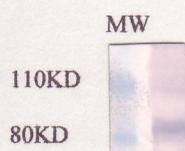

Western Blot analysis of OB0016 Chicken Anti-mouse HSL Abs (1:2000) on mouse adipose protein lysate (50ug).