

Product: Chicken Anti-Resistin, mouse and rat specific, affinity purified.

Catalog Number: OB0037

Lot Number: 80702

Size: 0.1mg in 0.1ml PBS with 0.02% Na Azide

Immunogen: Mouse resistin peptide VSSRGRGLASCPEGT.


Specification: OB0037 specifically binds to the mouse and rat resistin. The product does not cross react with human resistin.

Suggested Usage: The product can be used at 1:2,000 to 1:5,000 for immunoblot or ELISA testings. User needs to optimize the concentrations in assays. The product can be used in combination with GeneTel's anti-IgY conjugates such as RCYHRP (anti-IgY HRP) and RCYAP (anti-IgY AP) Western Blot analysis.

Short Term Storage: 2-8°C for one week.

Long Term Storage: -20°C as an undiluted liquid for up to two year.

Please feel free to call us if you have any question regarding to this product or any other GeneTel's product.

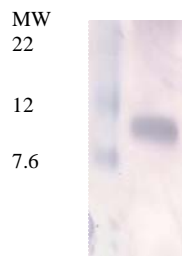

MW Ms Resistin

Western Blot analysis of OB0037 chicken anti-mouse resistin (1:2000) shows that the Abs binds strongly to the recombinant resistin protein (20ng) from PeproTech.