

Product: Chicken Anti-Resistin, human specific, affinity purified.

Catalog Number: OB0038

Lot Number: 80703

Size: 0.1mg in 0.1ml PBS with 0.02% Na Azide


Immunogen: Human resistin peptide AGMDNTGARCCRVQP.

Specification: OB0038 specifically binds to human resistin. The product may cross-react with bovine resistin but does not cross react with mouse or rat resistin.

Suggested Usage: The product can be used at 1:2,000 to 1:5,000 for immunoblot or ELISA tests. User needs to optimize the concentrations in assays. The product can be used in combination with GeneTel's anti-IgY conjugates such as RCYHRP (anti-IgY HRP) and RCYAP (anti-IgY AP) Western Blot analysis.

Short Term Storage: 2-8°C for one week.

Long Term Storage: -20°C as an undiluted liquid for up to two year.
Please feel free to call us if you have any question regarding to this product or any other GeneTel's product.

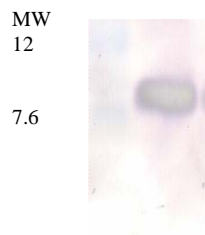

MW Hu resistin

Western Blot analysis of OB0038 chicken anti-human resistin (1:2000) shows that the Abs binds strongly to the human recombinant resistin protein (20ng) from PeproTech.