

GeneTel Laboratories, LLC *Leader in Proteomics Abs and Abs Arrays*
1202 Ann Street, Madison, WI 53713 Phone 877-248-4316 Fax 608-257-1385 Email info@genetel-lab.com

Product: Biotinylated Chicken Anti-Relaxin H2 B chain from human, affinity purified.

Catalog Number: OB0043

Lot Number: 130206

Size: 100ul at 0.5mg/ml in PBS with 0.02% Na Azide

Immunogen: Human Relaxin H2 B chain peptide HVGCTKRSLARFC.


Specification: OB0043 specifically binds to human Relaxin H2 B chain. The product may cross-react with relaxin H2 of other species.

Suggested Usage: The product can be used at 1:2,000 to 1:10,000 for ELISA tests in combination with OB0041 anti-H2 A chain. User needs to optimize the concentrations in assays. The product can be used in combination with GeneTel's streptavidin conjugates such as SAHRP (Streptavidin HRP) and SAAP (Streptavidin AP).

Short Term Storage: 2-8°C for one week.

Long Term Storage: -20°C as an undiluted liquid for up to two year.

Please feel free to call us if you have any question regarding to this product or any other GeneTel's product.

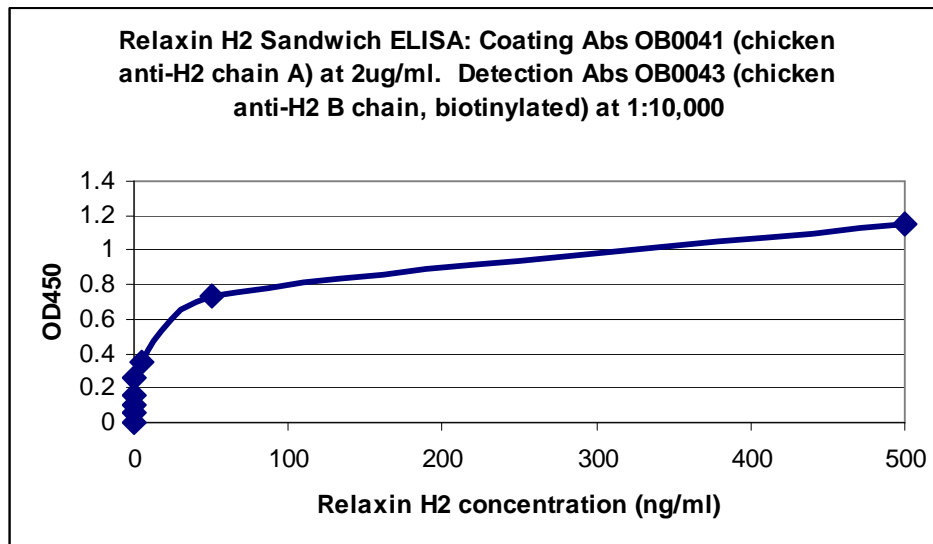

OB0043 in combination with OB0041 in ELISA testing for Relaxin H2