

GeneTel Laboratories, LLC *Leader in Proteomics Abs and Abs Arrays*
1202 Ann Street, Madison, WI 53713 Phone 877-248-4316 Fax 608-257-1385 Email info@genetel-lab.com

Product: Chicken Anti-H3 Relaxin A chain from human, affinity purified.

Catalog Number: OB0044

Lot Number: 130207

Size: 200ul at 0.5mg/ml in PBS with 0.02% Na Azide


Immunogen: Human Relaxin H3 A chain peptide KWGCSKSEISLCL.

Specification: OB0044 specifically binds to human H3 Relaxin A chain. The product may cross-react with H3 Relaxin of other species.

Suggested Usage: The product can be used at 1:2,000 to 1:5,000 for immunoblot or at 100-1:500 as coating Abs in H2 ELISA tests in combination with OB0046. User needs to optimize the concentrations in assays. The product can be used in combination with GeneTel's anti-IgY conjugates such as RCYHRP (anti-IgY HRP) and RCYAP (anti-IgY AP) Western Blot analysis.

Short Term Storage: 2-8°C for one week.

Long Term Storage: -20°C as an undiluted liquid for up to two year.
Please feel free to call us if you have any question regarding to this product or any other GeneTel's product.

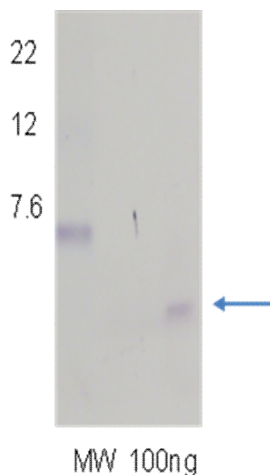

Western Blot analysis of OB0044 chicken anti-human H3 relaxin A chain (at 1:2000) shows that the Abs binds to the human recombinant relaxin H3 protein (100ng) of native form (non-denature on SDS-PAGE gel, about 5KDa) from PeproTech.